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Copyright, Eastman Kodak Company, 2004

Issued 3-04

KODAK Medical X-ray Film / General Purpose / Blue / 3530 / MXB

1) Description

KODAK Medical X-ray Film / General Purpose / Blue / 3530 / MXB is a full-speed, blue-sensitive medical x-ray screen film that is intended for general radiography. This film is coated on blue, 7-mil ESTAR Base with a base density of approximately 0.18. KODAK Medical X-ray Film / General Purpose / Blue features T-GRAIN emulsion technology that reduces the amount of screen-light crossover, resulting in improved image sharpness. It is designed for standard cycle machine processing.

2) Safelight

Use a KODAK GBX-2 Safelight Filter with a frosted 15-watt bulb or a KODAK LED Safelight located at least 4 feet from the film.

3) Storage and Handling

Handling -

Hands must be clean, dry and free of lotions, etc. Film should be handled carefully by the edges to avoid physical strains such as pressure, creasing, or buckling.

Storage -


Store unexposed film at 50 to 70°F (10 to 20°C), at 30 to 50 percent RH, and properly shielded from x-rays, gamma rays, or other penetrating radiation. Keep exposed film in a cool, dry place that is properly shielded from penetrating radiation. Process as soon as possible after exposure. Processed film should be stored at 60 to 80°F (16 to 27°C), at 30 to 50 percent RH.

4) Sensitometric Parameters

Relative Speed:	Measured at a density of 1.00 above gross fog.
Contrast:	Measured as slope of the line between densities of 0.25 and 2.00 above gross fog.
Gross Fog:	Density of film base plus processing fog.
Net Fog:	Processing fog.

5) Process Variations

Changes to speed, contrast, and fog as a result of temperature variation from normal are included in GRAPHS Section. Variations from bromide and iodide ions in KODAK RP X-OMAT Developer cause sensitometric speed effects that are significantly different for T-MAT Films than for conventional films.

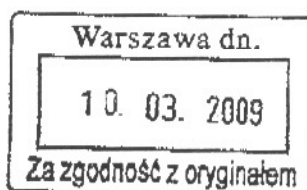AGFA Sp. z o.o.
Al. Jerozolimskie 195 A, 02-222 Warszawa
tel. 022 311-19-20, fax: 022 311-19-67
Regon 01000334, NIP 523-010-39-65

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
Data Sheet—5
T15019 Issued 3-04

Copyright, Eastman Kodak Company, 2004

T15019A 3-04
CHARACTERISTIC, For Publication

KODAK Medical X-ray Film / General Purpose / Blue / 3530 / M02
1/50 sec, Simulated Blue Screen, KODAK X-OMAT 5000 RA Processor,
Fresh KODAK RP X-OMAT Chemicals,
Diffuse Visual

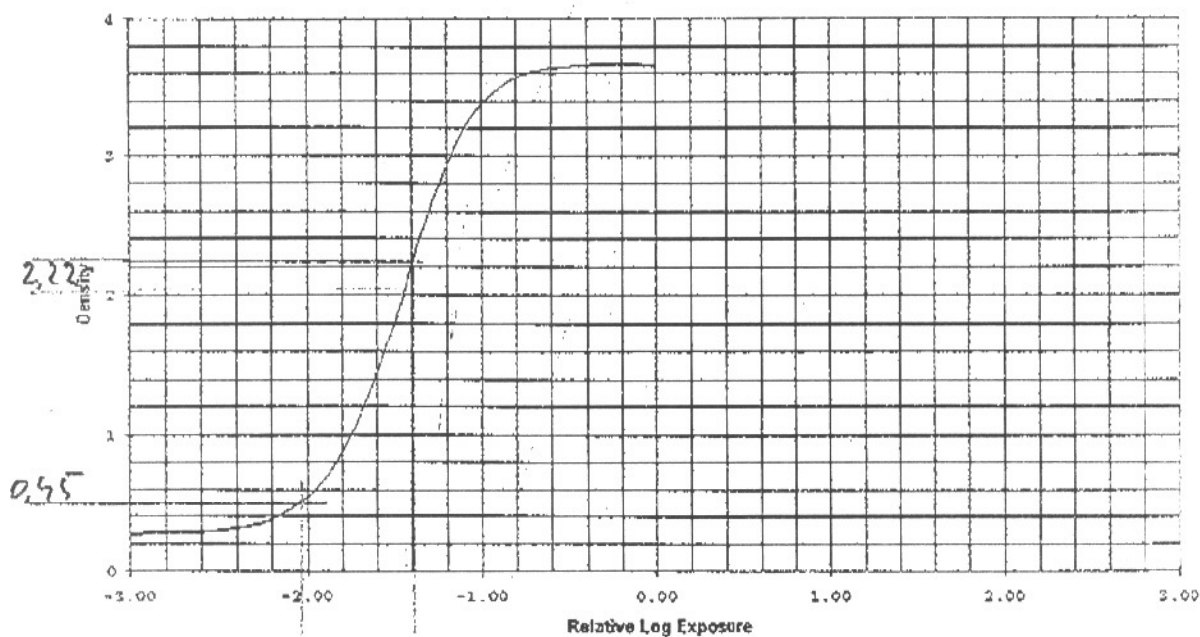Notice: While the data presented are typical of production coatings, they do not represent standards which must be met by Eastman Kodak Company. Varying storage, exposure, and processing conditions will affect results. The company reserves the right to change and improve the product characteristics at any time.

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